
Naipo

Instruction Manual

Shiatsu Pillow Massager with Heat for Back & Neck

Shenzhen Pinlong E-Commerce Co., Ltd.
Web: www.naipocare.com
E-mail: support@naipocare.com
Add: Huanan City Electronic Trading Building P05
Room 116, Longgang District, Shenzhen, China

Made in China

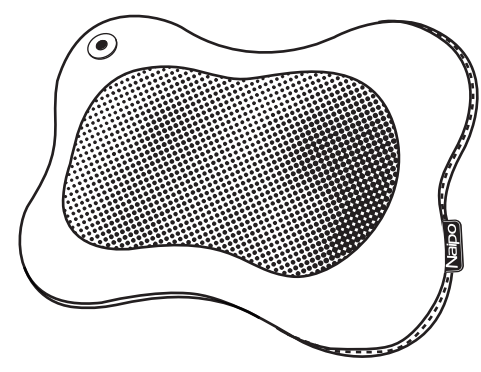Model: MGP-129M

Contents

English	01-04
Deutsch	05-08
Français	09-12
Español	13-16
Italiano	17-20

Contents	
English	01-04
Deutsch	05-08
Français	09-12
Español	13-16
Italiano	17-20

EN

This manual should be kept on hand for easy reference. It is your guide to a safe and efficient operation.

Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:


Danger

To reduce the risk of electric shock:
Always unplug this product from the electrical outlet immediately after use and before cleaning.


Warning

To reduce the risk of burns, fire, electric shock, or injury to persons:
1) The product should never be left unattended when plugged in. Unplug from outlet when it is not used and before putting on or taking off parts.
2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
3) Do not carry this product by its supply cord or use cord as a handle.
4) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
5) Do not drop or insert any object into any opening.
6) Do not use outdoors.
7) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
8) To disconnect, switch off all controls and remove plug from outlet.
9) Close supervision is necessary when this appliance is used by, or near children, invalids, or disabled persons to ensure that they do not misuse the product.
10) Keep the cord away from heated surfaces. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation.
11) The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
12) The appliance must only be supplied at safely extra low voltage corresponding to the marking on the appliance.
13) The appliance is only to be used with the power supply unit provided with the appliance.

Massage Position

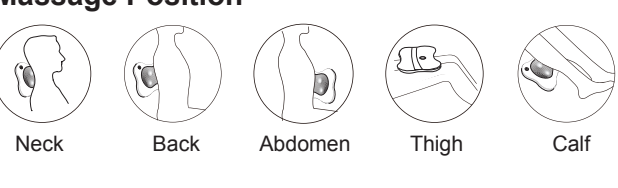Name and Function of The Machine Parts

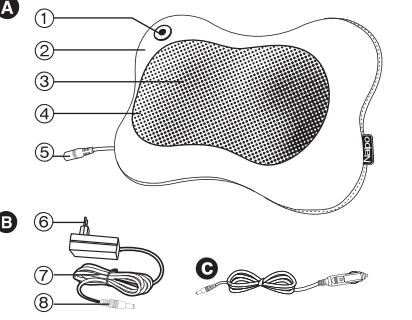- 1 Main Unit
- 2 Power Switch
- 3 Leather Cover
- 4 Massage Node
- 5 Lattice Cover
- 6 DC Connector
- 7 Power Adaptor
- 8 Plug
- 9 Cord
- 10 Appliance Plug
- 11 Vehicular Car Adaptor

How to Operate

1. Plug in the main plug of the car adaptor.
2. Press the Power Switch once. The product will start its kneading massage and heat therapy for 20 minutes.
3. Press the power switch again. This time, the product will stop its heat therapy while the kneading massage continues.
4. Press power switch and continue pressing for 3 seconds to stop all functions and massage session.

After use: When the product is not use, it should be unplugged.

To clean and maintain the product

- 1. Cleaning**
Before you clean, please make sure that the switch is in the off location, and the power plug has been pulled out from the socket. Please use dry cloth to clean the product.
- 2. Maintenance**
Keep the product clean and far from the wet environment.

Attention

Please do not use benzene or thinner and other aggressive detergents to clean the product or it may fade.

Contact Us

If you need any assistance or support, please contact us. We will do our best to provide you with a viable solution.
Web: www.naipocare.com
Email: support@naipocare.com

Technology parameter

Name: Massage Pillow
Model: MGP-129M
Rated Voltage: DC 12V = 2.5A
Rated Power: 30W

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