Model: MGS-N12CS

Shenzhen PinKing E-Commerce Co., Ltd. Web: www.naipocare.com Email: support@naipocare.com Add: Hunan City Electronic Trading Building P30 Room 114, Longgang District, Shenzhen, China. CE, RoHS, Made in China

Contents

English ----- 01-05
Deutsch ----- 06-10
Français ----- 11-15
Español ----- 16-20
Italiano ----- 21-25

Table with 2 columns: Model (MGS-N12CS) and specifications (Name, Input Voltage, Power, Default Time, Low Battery, Charging Time, Locking Time).

Frequently questions

Table with 3 columns: Question, Answer, Solution. Covers topics like broken, product not work, too slow and underpower, suddenly stop in working process, abnormal sound.


If the problem persists, please DO NOT USE continue and contact your distributor. Consult your doctor prior to use this product! You have a cardiac pacemaker, artificial joints or electronic implants.

Important safety warning Do not use the product with wet hands or put the wire hand. Do not use a power source other than the rated voltage. In case of electric shock, do not plug or unplug this product when your hands are wet.

Do not use the product if room temperature is higher than 40 degrees. Do not use the product in wet or too dusty environment to in case of electric shock or malfunction. Do not use in lightning weather and pull the power plug off the outlet.


Disassembly prohibited Do not disassemble the internal components of the product in case of product damage. In addition to our designated service technician, please do not make a decomposition, modification and maintenance management.

Cleaning and maintenance Maintenance Please use a mild detergent and a soft cloth to clean in daily maintenance. Do not use oil, organic solvents or other chemical agents and also a wire to scrub the blades.

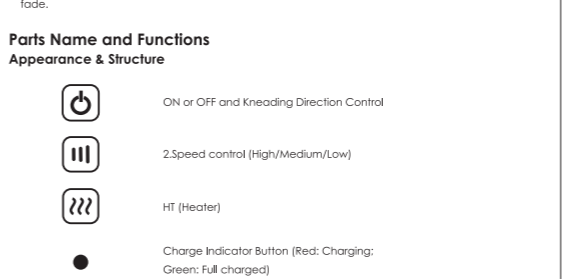Parts Name and Functions Control Buttons Description 1. Speed control 3 levels of speed, the default is high speed, press one time to low speed, press again middle speed, the high speed.

3. HI (Heater) You can choose to stop heating function when it is working, press one time, stop heating function, press again, heating function work again. Tips The product using the safety of the overheating protection device, it will stop when the temperature is overheating, you have to open again after cooling temperature.

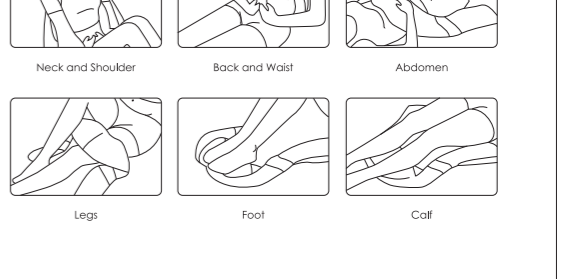Parts Name and Functions Control Buttons Description 1. Speed control 3 levels of speed, the default is high speed, press one time to low speed, press again middle speed, the high speed.

Table with 2 columns: Model (MGS-N12CS) and specifications (Name, Input Voltage, Power, Default Time, Low Battery, Charging Time, Locking Time, Nutritional Value).

Table with 2 columns: Model (MGS-N12CS) and specifications (Name, Input Voltage, Power, Default Time, Low Battery, Charging Time, Locking Time, Nutritional Value).

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